


Why Choose Servicebio Electronic Pipette Plus?

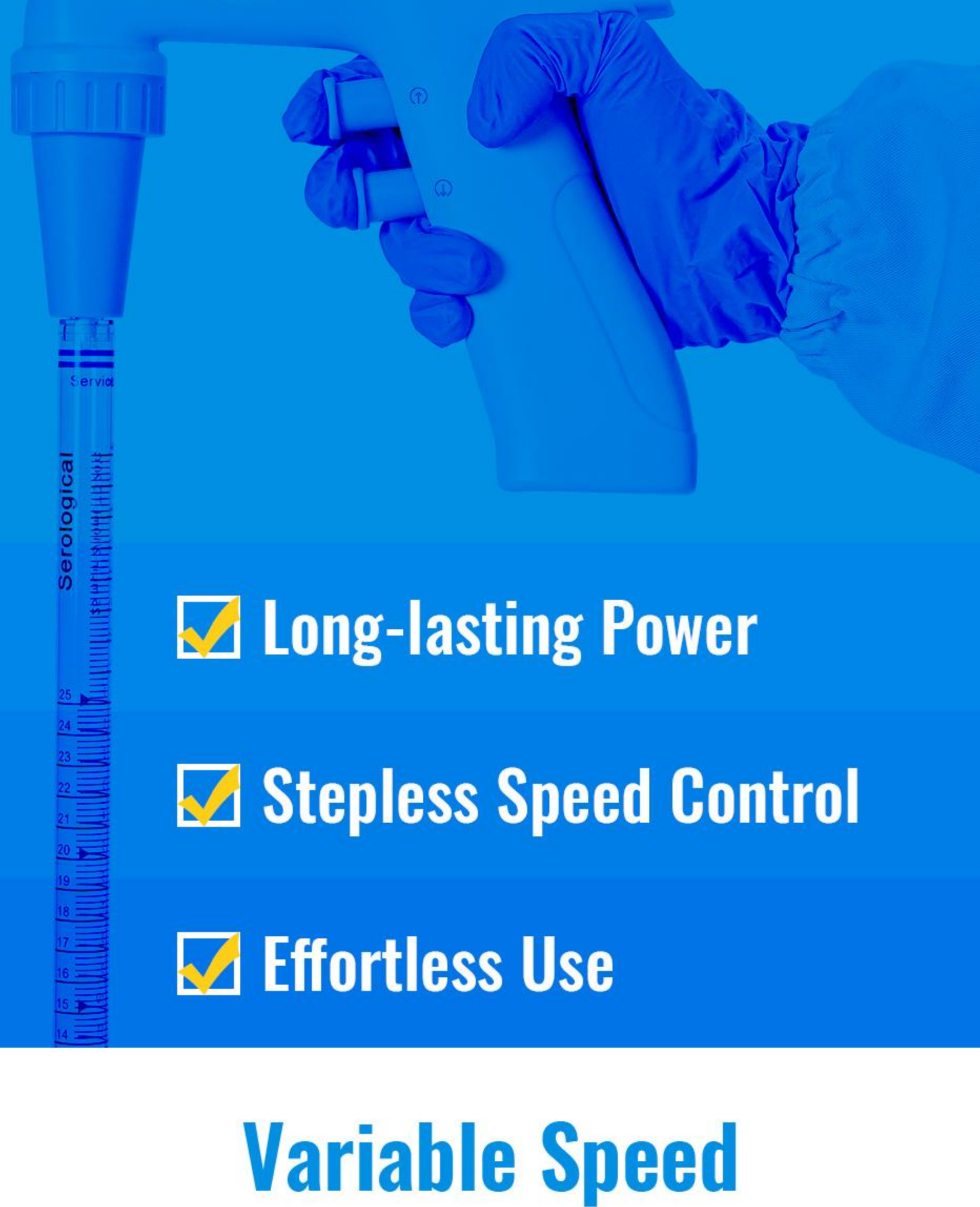

✓ Long-lasting Power

✓ Stepless Speed Control

✓ Effortless Use


Variable Speed

More Efficient Pipetting


Flexible speed adjustment by pressing force. No need to manually adjust the speed, easy and precise pipetting.



Gravity/Power Drainage

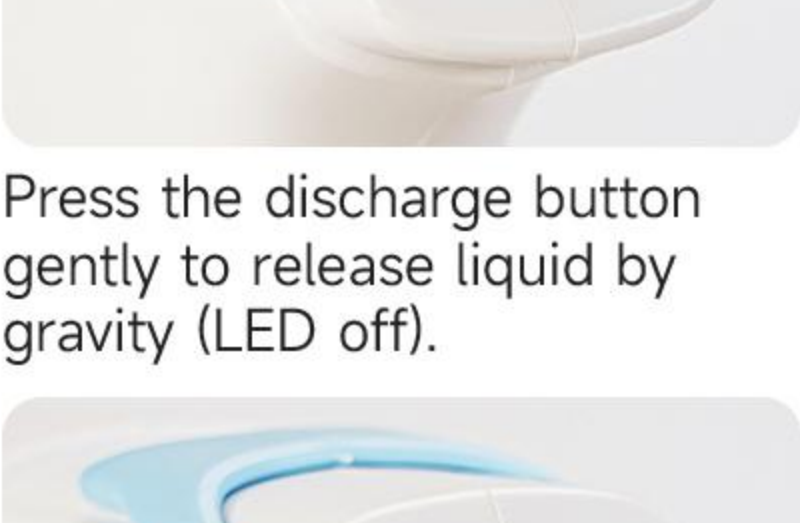

Press the discharge button gently to release liquid by gravity (LED off).

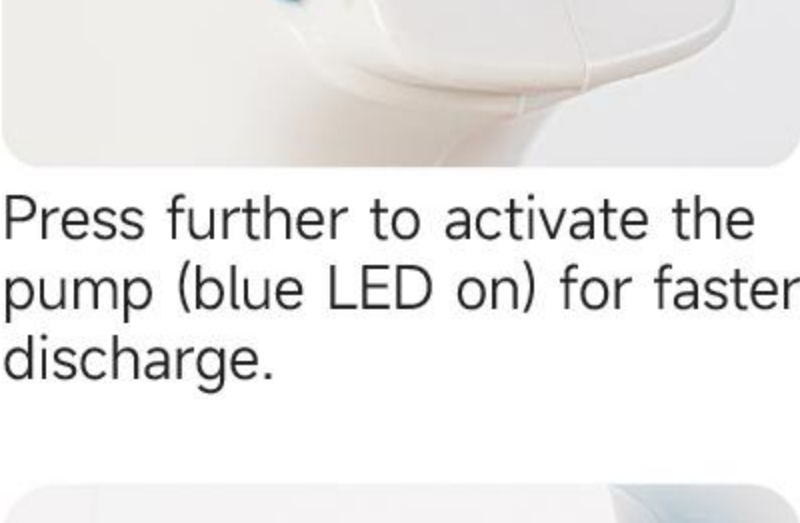

Press further to activate the pump (blue LED on) for faster discharge.










Aspiration

Drainage

Ultra-Long Battery life

1.6 h for Full Battery Charging

3 Weeks Intermittent Operation

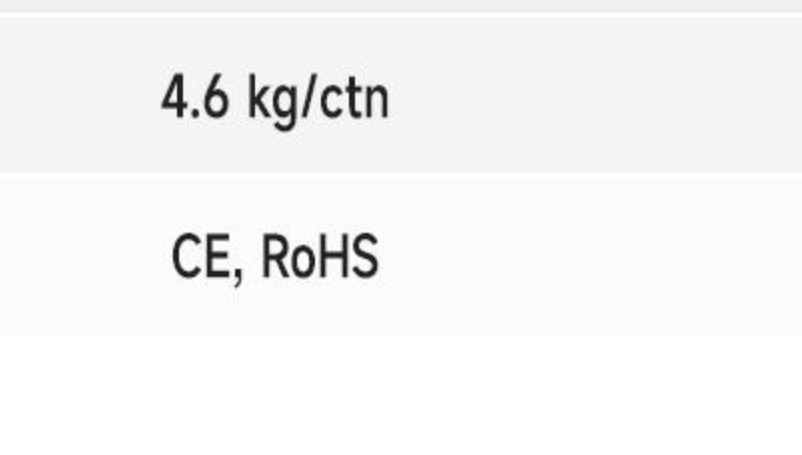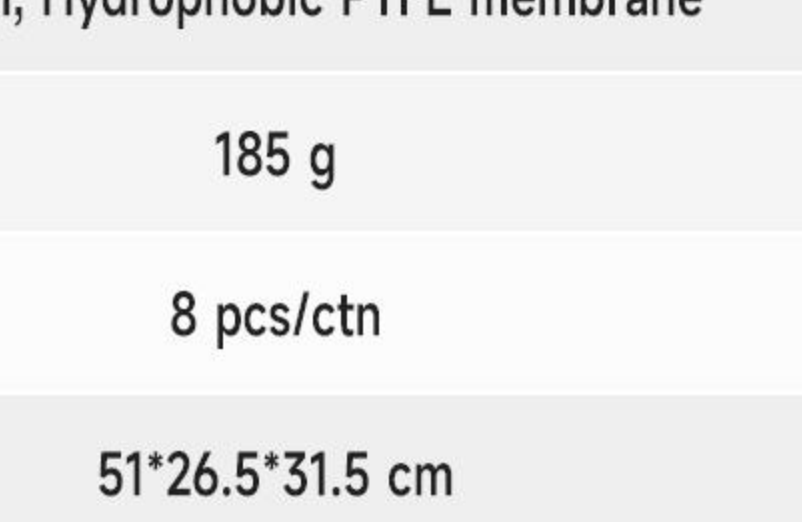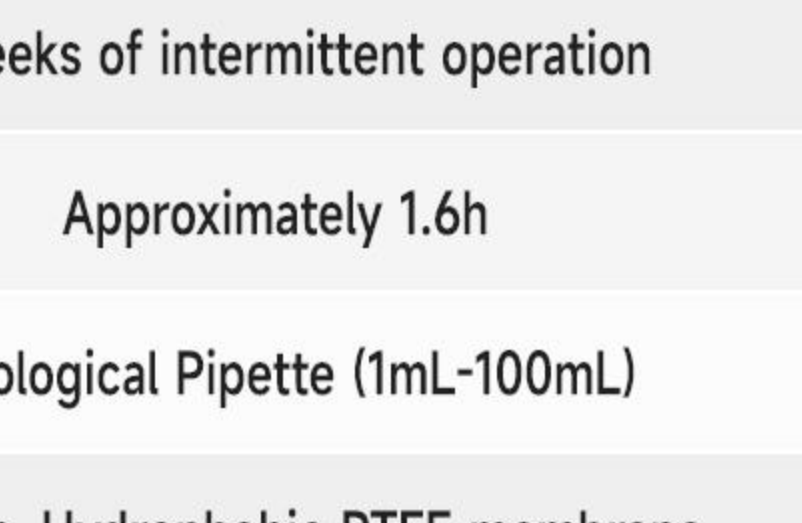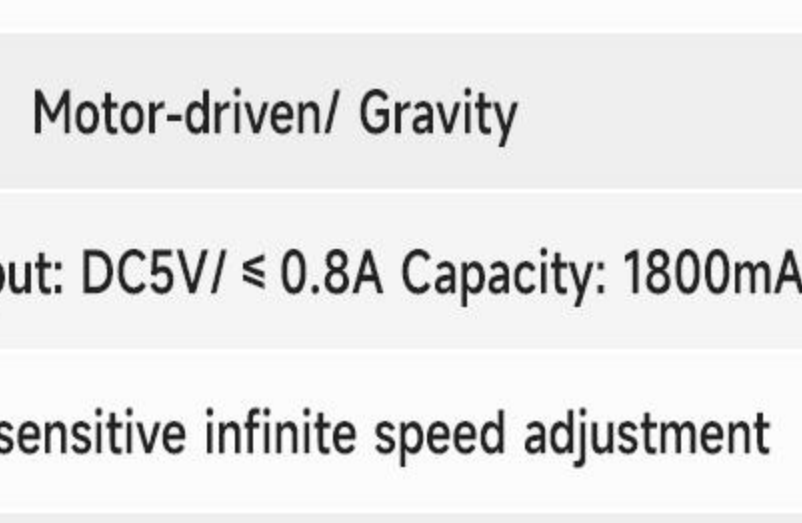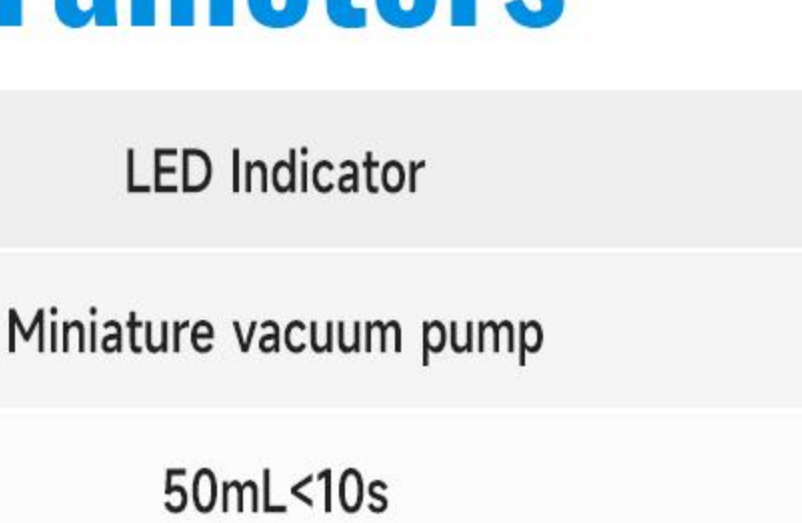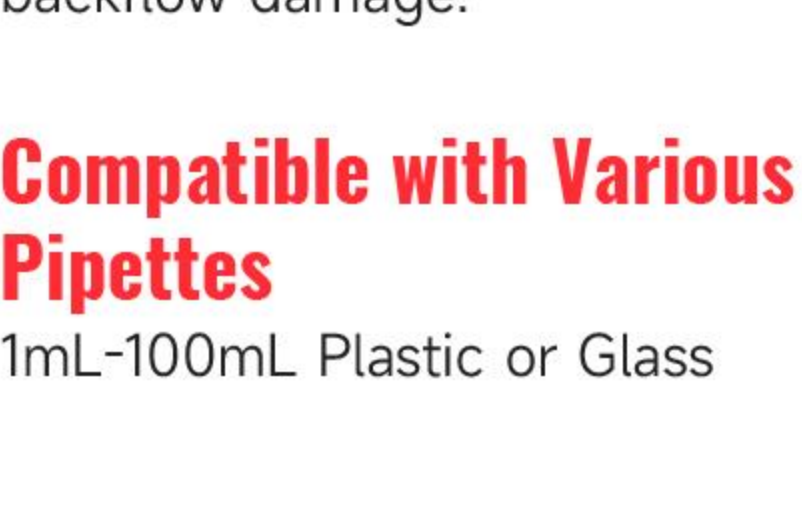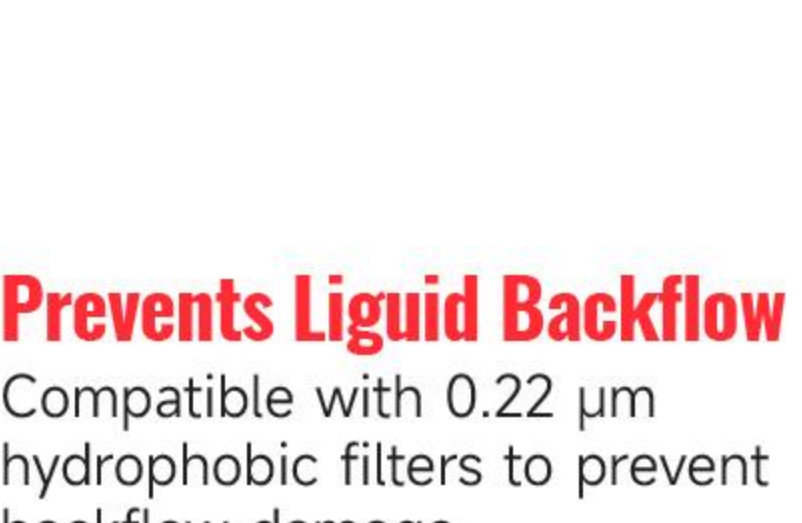

Product Package



Applications



Cell Culture Flasks



Centrifuge Tubes



Cell Culture Dish

Product Parameters

Battery Level Display	LED Indicator
Motor Type	Miniature vacuum pump
Aspiration Speed	50mL<10s
Drainage Speed	Motor-driven/ Gravity
Battery	Battery: Input: DC5V/ ≤0.8A Capacity: 1800mA
Speed Control	Pressure-sensitive infinite speed adjustment
Battery Life	3 weeks of intermittent operation
Charging Time	Approximately 1.6h
Compatible Pipette	Serological Pipette (1mL-100mL)
Filter	0.22µm, Hydrophobic PTFE membrane
Net Weight	185 g
Packaging	8 pcs/ctn
Carton Size (L*W*H)	51*26.5*31.5 cm
Gross Weight	4.6 kg/ctn
Certificate	CE, RoHS