

# Purite Fusion

Robust dual quality water purification system for the laboratories and life science sector



The Purite Fusion is a self-contained water purification unit that reliably delivers a steady supply of both Type I and Type II purified water from a mains supply to life science laboratories. Diagnostic uses include histology, cell and tissue culture, DNA sequencing and IVF.

## Additional features

- Dual water quality available 1-10MΩ.cm (Type II) and 18.2MΩ.cm (Type I)
- Optional Remote Dispense Pod
- Energy saving intelligent standby mode
- Dispense rate of up to 2 litres/min
- Water quality parameters, TOC, MΩ.cm, °C and flowrate displayed
- Selectable manual and volumetric dispense feature
- 5000m.wt cut off internal ultrafiltration for endotoxin, RNase/DNase removal
- Dual wavelength (185nm/254nm) UV irradiation
- Integral TOC indicator
- UV LED disinfection incorporated within the dispense head
- ECO option now available offering a significant reduction in water usage and waste
- WiFi enabled for remote monitoring and operation.

# Purite Fusion

Robust dual quality water purification system for the laboratories and life science sector


## Common Purite Features

Our Purite range of water purification systems are compact, robust, simple to use, easy to maintain and available in six standard models: Analyst, HP, Purewater 300 Package, Fusion, HPA 30 and Neptune Ultimate.

Common features of all our Purite systems include:

- Space-saving, dependable, bench top or wall mounted systems
- RO Removes > 98% minerals and > 99% bacteria
- Choice of production rates up to 48 l/hr
- Optional external storage tanks up to 100 litres
- RO Boost pump fitted as standard
- Installation kit and all consumables included for first few months operation
- LCD colour touch screen panel
- Visual and audible alarms included
- Utilises carbon pre-treatment, RO and deionisation
- USB port to download event data and upload software updates
- Integral 20 litre storage as standard (excludes Neptune Ultimate)
- Semi-automatic clean cycle.

## Purite Fusion Process Flow

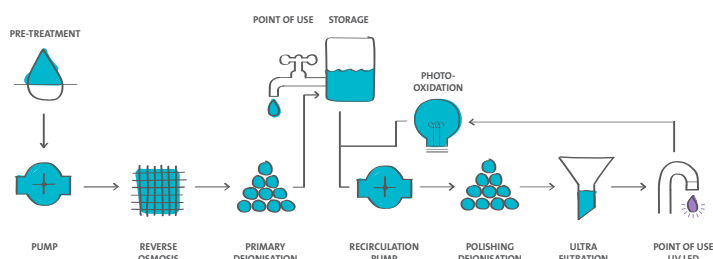

Visit online!

## PURITE

Veolia Water Purification Systems Ltd trading as Purite

Bandet Way Thame Oxfordshire OX9 3SJ United Kingdom

T +44 (0)1844 217141

E vtc.vwts.puritemail.all@veolia.com

W purite.com

## Technical Specifications

| Unit Specification                 | 160                                       | 320 |
|------------------------------------|-------------------------------------------|-----|
| Width (mm)                         | 440                                       |     |
| Depth (mm)                         | 560                                       |     |
| Height (mm)                        | 750                                       |     |
| Max shipping weight (kg)           | 36                                        | 41  |
| Max working weight (kg)            | 51                                        | 59  |
| Installation requirements          |                                           |     |
| Power                              | Single Phase, 110-230V, +/- 10%, 50/60 Hz |     |
| Feed water                         | Potable                                   |     |
| Maximum TDS (ppm)                  | 1000                                      |     |
| Minimum inlet pressure - psi (bar) | 10 (0.5)                                  |     |
| Maximum inlet pressure - psi (bar) | 30 (2.0)                                  |     |
| Feed water temperature             | 1-35°C                                    |     |
| Product outputs*                   |                                           |     |
| @ 10°C (l/hr)                      | 14.4                                      | 30  |
| @ 25°C (l/hr)                      | 24                                        | 48  |

\* Product outputs based on a feed water pressure of 4 bar

| System Specification          |                                                                                |
|-------------------------------|--------------------------------------------------------------------------------|
| Pure water storage            | 20 litre storage tank as standard<br>(External 50 & 100 litre tanks available) |
| Display panel                 | LCD – Colour touch screen                                                      |
| Pre-treatment cartridge       | ✓                                                                              |
| Reverse osmosis               | ✓                                                                              |
| Deionisation cartridge        | ✓                                                                              |
| Internal filtration           | Ultrafiltration                                                                |
| Point of use                  | UV / LED                                                                       |
| UV lamp                       | 185nm / 254nm                                                                  |
| Recirculation pump            | ✓                                                                              |
| Ultrapure polishing cartridge | ✓                                                                              |

| Treated Water Specification | High Purity Dispense      | Purified Water Storage Tank |
|-----------------------------|---------------------------|-----------------------------|
| Inorganics                  | up to 18.2MΩ.cm           | > 1MΩ.cm                    |
| pH*                         | Neutral                   |                             |
| Bacteria                    | < 0.1cfu/ml               | -                           |
| Organics – TOC (ppb)        | < 5                       | < 50                        |
| Particles                   | < 0.005µm                 | -                           |
| Endotoxins                  | < 0.001EU/ml              | -                           |
| DNases                      | < 4pg/µl                  | -                           |
| RNases                      | < 0.01ng/ml               | -                           |
| Dispense modes              | Manual - Latched - Volume |                             |
| Dispense flow rate          | up to 2.0 l/min           | -                           |

\* pH of stored water may decrease due to absorption of free carbon dioxide